

### Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell


Catalog No.	Size
CHEK-ATF167	$2 \times (1 \text{ vial contains } \sim 5 \times 10^6 \text{ cells})$

### • Description

The Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell was engineered to not only express STAT3 signaling response element, but also express the receptors full length human IL-22 R alpha 1 (Gene ID: 58985) and IL-10 R beta (Gene ID: 3588). When stimulated with human IL-22 protein, receptor-mediated signaling can drive STAT3-mediated luminescence. Neutralization of biological effect of human IL-22 protein by corresponding antibody results in a decrease in luminescence.

#### • Application

- Screen for neutralizing antibodies blocking the stimulation of human IL-22 protein.
- Bioactivity detection of human IL-22 fusion protein.

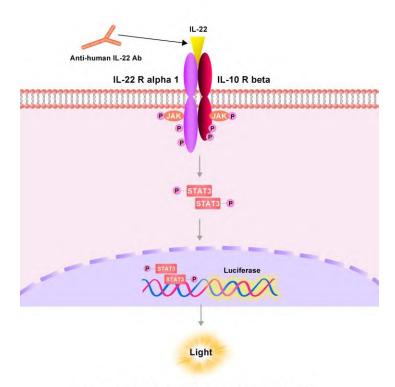

Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell



#### • Cell Line Profile

Cell line	Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell
Host Cell	HEK293
Property	Adherent
Complete Growth Medium	DMEM + 10% FBS
Selection Marker	Puromycin (2 μg/mL) + Hygromycin (40 μg/mL) + Zeocin (20 μg/mL)
Incubation	37°C with 5% CO <sub>2</sub>
Doubling Time	22-24 hours
Transduction Technique	Lentivirus

#### • Materials Required for Cell Culture

- DMEM medium (Gibco, Cat.No.11965-092)
- Fetal bovine serum (CellMax, Cat.No.SA211.02)
- Puromycin (InvivoGen, Cat.No.ant-pr-5b)
- Hygromycin B (Invitrogen, Cat.No.10687010)
- Zeocin (Invitrogen, Cat.No.R25001)
- 0.25% Trypsin-EDTA (1X), Phenol Red (Gibco, Cat.No.25200-056)
- Penicillin-Streptomycin (Gibco, Cat.No.15140-122)
- Phosphate Buffered Saline (1X) (HyClone, Cat.No.SH30256.01)
- Complete Growth Medium: DMEM + 10% FBS, 1% P/S
- Culture Medium: DMEM + 10% FBS, Puromycin (2 μg/mL), Hygromycin (40 μg/mL), Zeocin (20 μg/mL),
  1% PS
- Freeze Medium: 90% FBS, 10% (V/V) DMSO
- T-75 Culture flask (Corning, 430641)
- Cryogenic storage vials (SARSTEDT, 72.379.007)
- Thermostat water bath
- Centrifuge
- Luna cell counter (Logos Biosystems, LUNA-II)
- CO<sub>2</sub> Incubator (Thermo, 3111)
- Biological Safety Cabinet (Thermo, 1389)



#### • Recovery

- 1. Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the cap out of the water. Thawing should be rapid (approximately 2 minutes).
- 2. Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by spraying with 70% ethanol. All the operations from this point on should be carried out under strict aseptic conditions.
- 3. Transfer the vial contents to a centrifuge tube containing 4.0 mL complete growth medium and spin at approximately 1000 rpm for 5 minutes.
- 4. Resuspend cell pellet with 5 mL complete growth medium and transfer the cell suspension into T-75 flask containing 10-15 mL of pre-warmed complete growth medium.
- 5. Incubate at 37°C with 5% CO<sub>2</sub> incubator until the cells are ready to be split.

#### • Subculture

- 1. Remove and discard culture medium.
- 2. Wash the cells once with sterile PBS.
- 3. Add 2 mL of 0.25% trypsin to cell culture flask. Place the flask at 37°C for 2-3 minutes, until 90% of the cells have detached.
- 4. Add 6.0 to 8.0 mL of culture medium and aspirate cells by gently pipetting.
- 5. Add appropriate aliquots of the cell suspension to new culture vessel.
- 6. Incubate at 37°C with 5% CO<sub>2</sub> incubator.

**Subcultivation Ratio:** A subcultivation ratio of 1:6 to 1:10 is recommended.

**Medium Renewal:** Every 2 to 3 days.

**Note:** After recovery for 1-2 generations with the complete growth medium not containing the selection marker, if the cell state is well, changing to the culture medium containing the selection marker.



#### • Cryopreservation


- 1. Remove and discard spent medium.
- 2. Detach cells from the cell culture flasks with 0.25% trypsin.
- 3. Centrifuge at 1000 rpm for 5 min at RT to pellet cells.
- 4. Resuspend the cell pellets with complete growth medium and count viable cells.
- 5. Centrifuge at 1000 rpm for 5 min at RT and resuspend cells in freezing medium to a concentration of  $5 \times 10^6$  to  $1 \times 10^7$  cells/mL.
- 6. Aliquot into cryogenic storage vials. Place vials in a programmable cooler or an insulated box placed in a -80°C freezer overnight, then transferring to liquid nitrogen storage.

#### • Storage

- Product format: Frozen
- Storage conditions: Liquid nitrogen immediately upon receipt



### • Receptor Assay

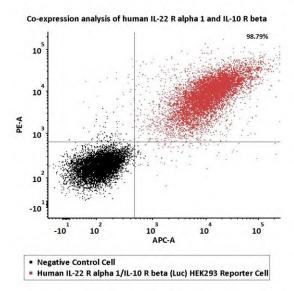

**Fig1.** Co-expression analysis of human IL-22 R alpha 1 and IL-10 R beta on Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell by FACS. Cell surface staining was performed on Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell or negative control cell using anti-human IL-22 antibody followed by staining with PE anti-human IgG Fc antibody and APC-labeled anti-IL-10 R beta antibody.

#### • Signaling Bioassay


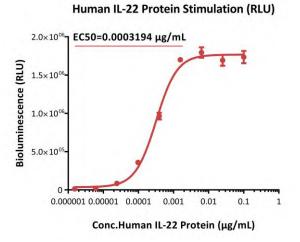


Fig2. Response to human IL-22 protein (RLU). The Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell was stimulated with serial dilutions of human IL-22 protein (Cat.No.IL2-H524a). The EC50 was approximately  $0.0003194 \, \mu g/mL$ .



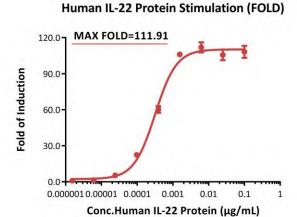

**Fig3. Response to human IL-22 protein (FOLD).** The Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell was stimulated with serial dilutions of human IL-22 protein (Cat.No.IL2-H524a). The max induction fold was approximately 111.91.

### Application


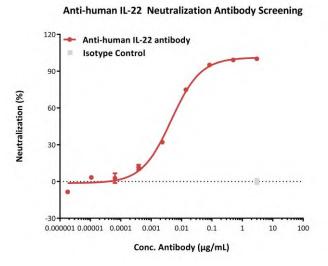


Fig4. Inhibition of human IL-22 protein-induced reporter activity by anti-human IL-22 neutralizing antibody. The Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell was incubated with serial dilutions of antibodies in the presence of human IL-22 protein (Cat.No.IL2-H524a) with a final concentration of  $0.0005~\mu\text{g/mL}$ . The EC50 of anti-human IL-22 neutralizing antibody is approximately  $0.004567~\mu\text{g/mL}$ .




#### • Passage Stability

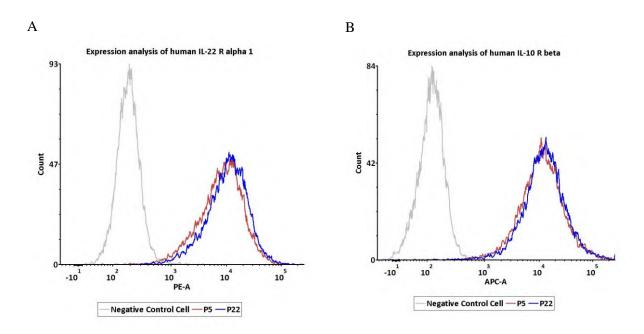

Passage	MFI for IL-22 R alpha 1 (PE)	MFI for IL-10 R beta (APC)
P5	8368.39	10859.46
P22	10555.09	12303.22

**Fig5. Passage stability analysis of receptors expression by FACS.** Flow cytometry surface staining of human IL-22 R alpha 1 and IL-10 R beta on Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell demonstrates consistent mean fluorescent intensity across across passage 5-22. (A) Human IL-22 R alpha 1 expression analysis. (B) Human IL-10 R beta expression analysis.



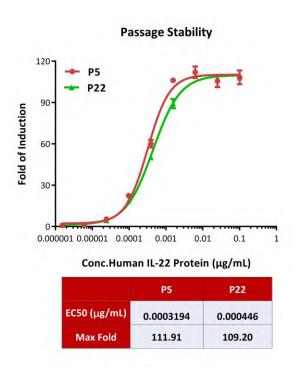

**Fig6. Passage stability analysis by Signaling Bioassay.** The continuously growing Human IL-22 R alpha 1/IL-10 R beta (Luc) HEK293 Reporter Cell was stimulated with serial dilutions of IL-22 protein. IL-22 protein stimulated response demonstrates passage stabilization (fold induction and EC50) across passage 5-22.

#### • License Disclosure

This reporter cell is provided for research use only. This license does not permit you to share, distribute, sell, sublicense, or otherwise make this reporter cell available for use to other laboratories, departments, research institutions, hospitals, universities, or biotech companies. The license does not permit modification of this reporter cell in any way. Inappropriate use or distribution of this reporter cell will result in revocation of the license. Modifications of this cell line, transfer to another facility, or commercial use of the cells may require a separate license and additional fees. AcroBiosystems does not warrant the suitability of this reporter cell for any particular use, and does not accept any liability in connection with the handling or use of this reporter cell.

#### • Related Products

<u>Products</u> Human IL-22 Protein, His Tag Cat.No.

IL2-H524a