Catalog # LDR-H82E7

ACCO

Synonym

LDLR,FH,FHC,LDLCQ2

Source

Biotinylated Human LDL R, His, Avitag(LDR-H82E7) is expressed from human 293 cells (HEK293). It contains AA Ala 22 - Arg 788 (Accession # <u>P01130-1</u>). Predicted N-terminus: Ala 22

Molecular Characterization

LDL R(Ala 22 - Arg 788) P01130-1 Poly-his Avi

This protein carries a polyhistidine tag at the C-terminus, followed by an Avi tag (AvitagTM).

The protein has a calculated MW of 88.5 kDa. The protein migrates as 110-115 kDa and 130-150 kDa when calibrated against <u>Star Ribbon Pre-stained Protein</u> <u>Marker</u> under reducing (R) condition (SDS-PAGE) due to glycosylation.

Labeling

Biotinylation of this product is performed using Avitag[™] technology. Briefly, the single lysine residue in the Avitag is enzymatically labeled with biotin.


Protein Ratio

Passed as determined by the HABA assay / binding ELISA.

Endotoxin

Less than 1.0 EU per μg by the LAL method.

SDS-PAGE

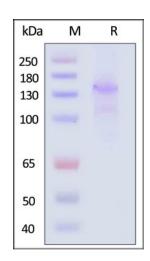

Biotinylated Human LDL R, His, Avitag on SDS-PAGE under reducing (R) condition. The gel was stained with Coomassie Blue. The purity of the protein

Purity

>80% as determined by SDS-PAGE.

Formulation

Lyophilized from 0.22 μ m filtered solution in PBS, pH7.4 with trehalose as protectant.

Contact us for customized product form or formulation.

Reconstitution

Please see Certificate of Analysis for specific instructions.

For best performance, we strongly recommend you to follow the reconstitution protocol provided in the CoA.

Storage

For long term storage, the product should be stored at lyophilized state at -20°C or lower.

Please avoid repeated freeze-thaw cycles.

This product is stable after storage at:

- -20°C to -70°C for 12 months in lyophilized state;
- -70°C for 3 months under sterile conditions after reconstitution.

is greater than 80% (With Star Ribbon Pre-stained Protein Marker).

Bioactivity-ELISA






Biotinylated Human LDL R Protein, His,Avitag™



Catalog # LDR-H82E7

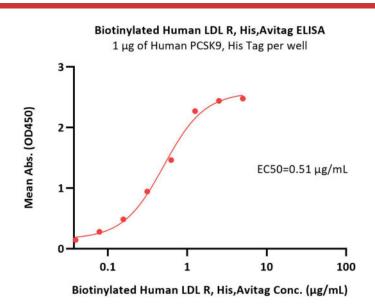

Immobilized Human PCSK9, His Tag (Cat. No. PC9-H5223) at 10 μ g/mL (100 μ L/well) can bind Biotinylated Human LDL R, His,Avitag (Cat. No. LDR-H82E7) with a linear range of 0.02-0.625 μ g/mL (QC tested).

Background

Low-Density Lipoprotein (LDL) Receptor is also known as LDLR, FH, FHC, LDLCQ2, and is a mosaic protein of ~840 amino acids (after removal of signal peptide) that mediates the endocytosis of cholesterol-rich LDL. It is a cell-surface receptor that recognizes the apoprotein B100 which is embedded in the phospholipid outer layer of LDL particles. The receptor also recognizes the apoE protein found in chylomicron remnants and VLDL remnants (IDL). It belongs to the Low density lipoprotein receptor gene family. LDL receptor complexes are present in clathrin-coated pits (or buds) on the cell surface, which when bound to LDL-cholesterol via adaptin, are pinched off to form clathrin-coated vesicles inside the cell. This allows LDL-cholesterol to be bound and internalized in a process known as endocytosis and prevents the LDL just diffusing around the membrane surface. This occurs in all nucleated cells (not erythrocytes), but mainly in the liver which removes ~70% of LDL from the circulation. Synthesis of receptors in the cell is regulated by the level of free intracellular cholesterol; if it is in excess for the needs of the cell then the transcription of the receptor gene will be inhibited. LDL receptors are translated by ribosomes on the endoplasmic reticulum and are modified by the Golgi apparatus before travelling in vesicles to the cell surface. LDL is directly involved in the development of atherosclerosis, due to accumulation of LDL-cholesterol in the blood. Atherosclerosis is the process responsible for the majority of cardiovascular diseases.

Clinical and Translational Updates

Please contact us via <u>TechSupport@acrobiosystems.com</u> if you have any question on this product.



